



**THE No. 1 System
for Solvent Safety**

Arnold Ball-Valve Funnel V2.0 with Hinged Lid

Strong Protection when collecting
Liquid Waste!



Clever Technology, Strong Protection!

Arnold Safety Funnel

✓ Improved Closing

Comfortable handling: easy, quick opening and closing with "click"-closure.

✓ Hinged Lid with Sealing Ring

For secure closure, when the funnel is not in use.

✓ Stable Hinge

A double-screw hinge, constructed of rust-free stainless steel: for long life and stability.

✓ Optimized Spray Protection

Air ducts directed downwards ensure for secure ventilation and cause any spray created to be returned to the funnel.

✓ Optimized Dirt Sieve

For collecting mixing magnets or large dirt particles. Laboratory bottles can also be placed here to dry off. The hole size has been optimized, in order to improve the running-off of chemicals.

Both sieve and spray protection can be removed for cleaning.

✓ Ball-Valve

The valve only opens up during filling, and then closes automatically afterwards - for a big "Plus" regarding safety!

✓ Various Screw-Threads

Suitable for different sizes of collecting canister (see Table on Page 3).

✓ Safety Lance

Ensures for safe and good flowing of sensitive chemicals, without any "stuttering".

✓ Robust Design

The new design provides for a high degree of stability as well as long life, and prevents any deformation of the funnel from occurring.

Information for Ordering

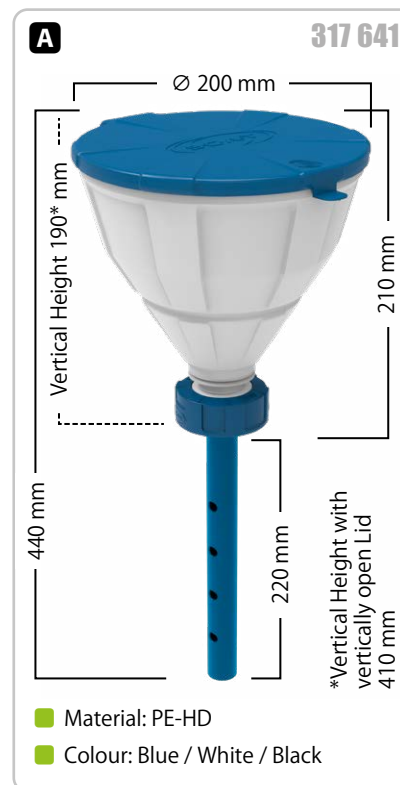


Fig.	Part No.	Screw-Thread of Container	Lance	Electrically conductive	Delivery with Adapter
-	317 642	GL 45	●	-	-
-	317 622	GL 45	●	●	-
-	317 638	GL 45	-	●	-
-	317 649	S 50	●	-	●
-	317 629	S 50	●	●	●
-	317 644	S 51	●	-	●
-	317 624	S 51	●	●	●
-	317 632	B 53	●	-	●
-	317 645	S 55	●	-	-
-	317 625	S 55	●	●	-
A	317 641	S 60 / 61	●	-	-
B	317 621	S 60 / 61	●	●	-
-	317 646	S 65	●	-	●
-	317 626	S 65	●	●	●
-	317 648	S 70 / 71	●	-	●
-	317 628	S 70 / 71	●	●	●
-	317 647	B 83	●	-	●
-	317 627	B 83	●	●	●
-	317 643	S 90	●	-	●
-	317 623	S 90	●	●	●
-	317 635	S 95	●	●	●
-	317 630	Dual Screw-Thread R2" BSP/ G2" + 2" TriSure	-	●	-
Models without ball valve, for highly viscous liquids					
-	317 633	GL 45	-	●	-
-	317 634	S 60 / 61	-	●	-



Information about "Delivery with Adapter":

The safety funnel "Arnold" is available in various container screw-thread sizes. In order that we may produce in an environmentally-friendly and economic manner, certain variants are delivered together with a screw-thread adapter. Those to which this applies can be seen in the column "Delivery with Adapter".



* The total height can vary, dependent upon the thread adapter with which the funnel has been delivered.





SCAT - Europe GmbH

Waldecker Straße 7
64546 Mörfelden-Walldorf
Germany



+49 (0) 61 05 / 30 55 86 - 0



+49 (0) 61 05 / 30 55 86 - 99



info@scat-europe.com



www.scat-europe.com



SCAT - Americas, Inc.

1704 East Blvd., Suite 101
Charlotte, NC 28203
USA



+1.704.438.9477



+1.704.373.2603



info@scat-americas.com



www.scat-americas.com

